


Record Temperature Readings (Manual)

1. Click on 'Storage Units' from the 'Clinic Tools' dropdown in the main menu.
2. Click on the 'Log' dropdown.

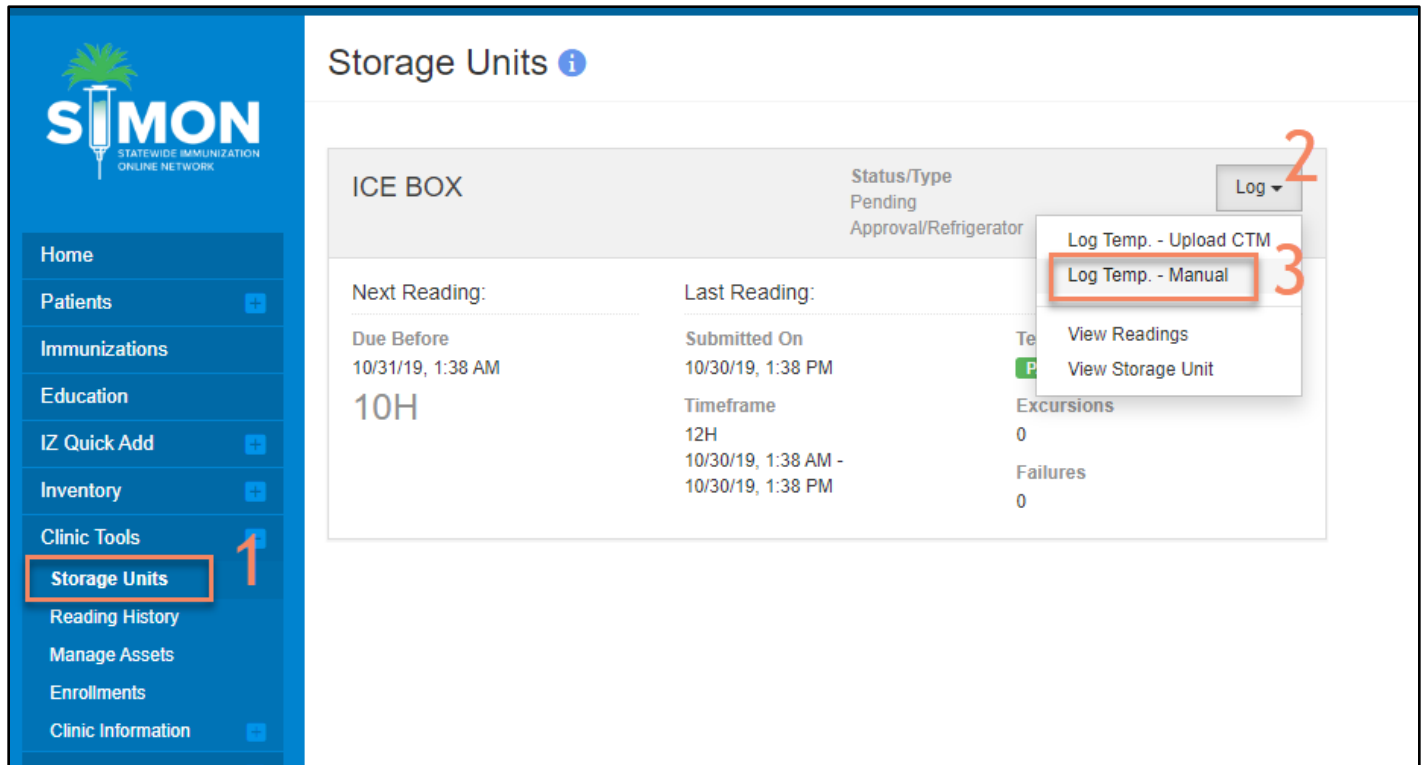

The screenshot displays the SIMON interface for 'Storage Units'. On the left is a blue navigation menu with 'Storage Units' highlighted under 'Clinic Tools'. The main content area shows details for an 'ICE BOX' with a status of 'Pending Approval/Refrigerator'. A 'Log' dropdown menu is open, showing options: 'Log Temp. - Upload CTM', 'Log Temp. - Manual' (highlighted with a red box and number 3), 'View Readings', and 'View Storage Unit'. The 'Log' dropdown is also highlighted with a red box and number 2. The 'Log Temp. - Manual' option is highlighted with a red box and number 3. The 'Storage Units' menu item is highlighted with a red box and number 1.

| ICE BOX | | Status/Type | Log |
|-------------------|---------------------|--------------------------------------|-------------------|
| | | Pending Approval/Refrigerator | Log |
| Next Reading: | Last Reading: | | |
| Due Before | Submitted On | Te | View Readings |
| 10/31/19, 1:38 AM | 10/30/19, 1:38 PM | P | View Storage Unit |
| 10H | Timeframe | Excursions | |
| | 12H | 0 | |
| | 10/30/19, 1:38 AM - | Failures | |
| | 10/30/19, 1:38 PM | 0 | |

3. Select 'Log Temp. – Manual'.

4. Enter the following information: Date, Time, Temperature, Min-Max Temperatures.

Log Manual Temperature Reading ?

| | | |
|--|---|---------------------------------|
| Storage Units / ICE BOX (Refrigerator) | Thermometer (Type) THERMOMETER (CTM) | Unit Status Pending Approval |
|--|---|---------------------------------|

The associated storage unit is not approved. Vaccines cannot be stored in the unit until approved.

Temperature must be logged as Fahrenheit.

| | | |
|---|---|---|
| Date * <input type="text" value="10/30/2019"/> | Time * <input type="text" value="HH:MM AM/PM"/> | Temperature * <input type="text" value="TEMPERATURE"/> °F |
| Min Temperature * <input type="text" value="MIN"/> °F | Max Temperature * <input type="text" value="MAX"/> °F | |

5. Click 'Log Temperature'.

| | |
|--|--|
| Next Reading 10/31/19, 1:38 AM | Clinic Notes Expand <input type="button" value="+ Add"/> There are currently no notes entered for this clinic |
|--|--|


Storage Unit

[Edit Storage Unit](#)

[Temperature Readings](#)

Log Temperature

6. A green success message is displayed confirming your changes were saved.

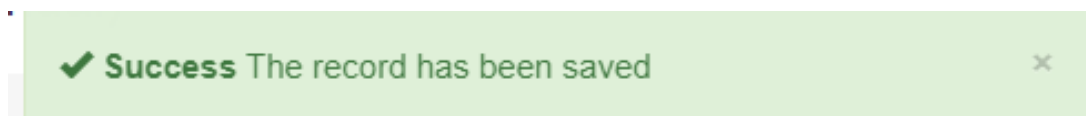

7. If the reading triggers a failure, you will receive an error message.

- Once you have added a reading successfully, you are redirected to the 'Storage Units' screen and the most recently added readings are displayed.

Temperature Readings i

Storage Units / ICE BOX (Refrigerator)

Thermometer (Type)
THERMOMETER (CTM)

Unit Status
Pending Approval

Showing 1 to 2 of 2 entries

| Submitted On | Reading Type | Timeframe | Failures | Min | Max | Action |
|--------------------|--------------|-----------|----------|--------|--------|----------------------|
| 10/30/19 - 2:46 PM | MANUAL | 12H | 0 | 35.0°F | 41.0°F | VIEW |
| 10/30/19 - 1:38 PM | MANUAL | 12H | 0 | 35.0°F | 41.0°F | VIEW |

Showing 1 to 2 of 2 entries

← Previous 1 Next →

- Click on 'View' to look at reading details.

Temperature Reading - 10/14/19, 12:31 PM (Manual) i

Storage Units / ICE BOX (Refrigerator)

Thermometer (Type)
THERMOMETER 1 (CTM)

Unit Status
Active

Next Reading
10/15/19, 12:31 AM

Summary

| | | |
|--|-----------------------|---------------------------------|
| Date 10/14/19, 12:31 PM | Temperature 40.0°F | Min/Max Temp 36.0°F / 41.0°F |
| Total Timeframe Total Duration: 12H From: 10/14/19, 12:31 AM To: 10/14/19, 12:31 PM | Excursions 0 | Failures 0 |

[Action ▼](#)

Excursions

Showing 0 to 0 of 0 entries

| Date | Temperature |
|-------------------|-------------|
| NO DATA AVAILABLE | |

Showing 0 to 0 of 0 entries

Storage Unit

- [Edit Storage Unit](#)
- [Temperature Readings](#)
- [Log Temperature](#)

Associated Thermometer

- [Edit Thermometer](#)
- [Calibrate Thermometer](#)

← Previous Next →